




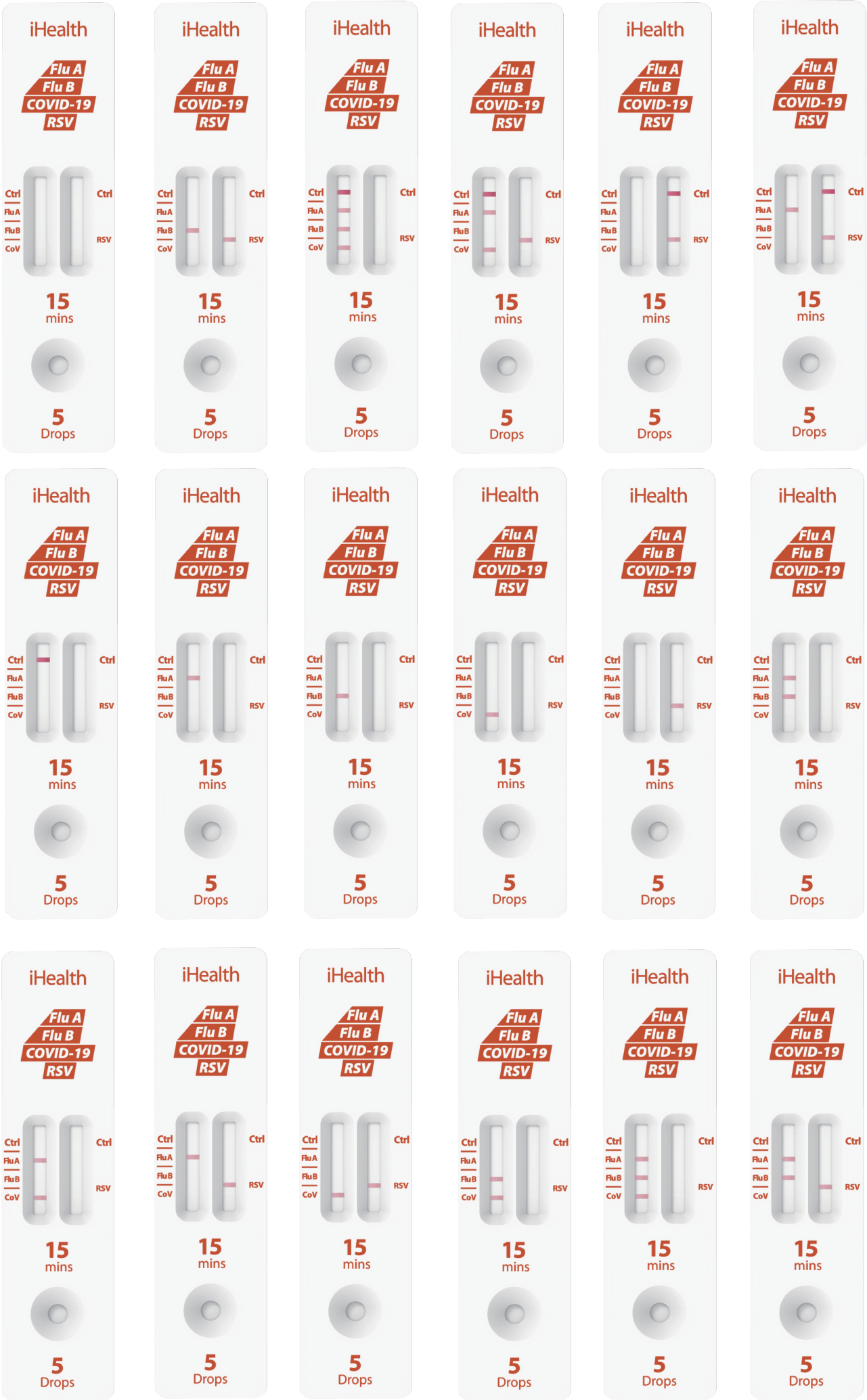

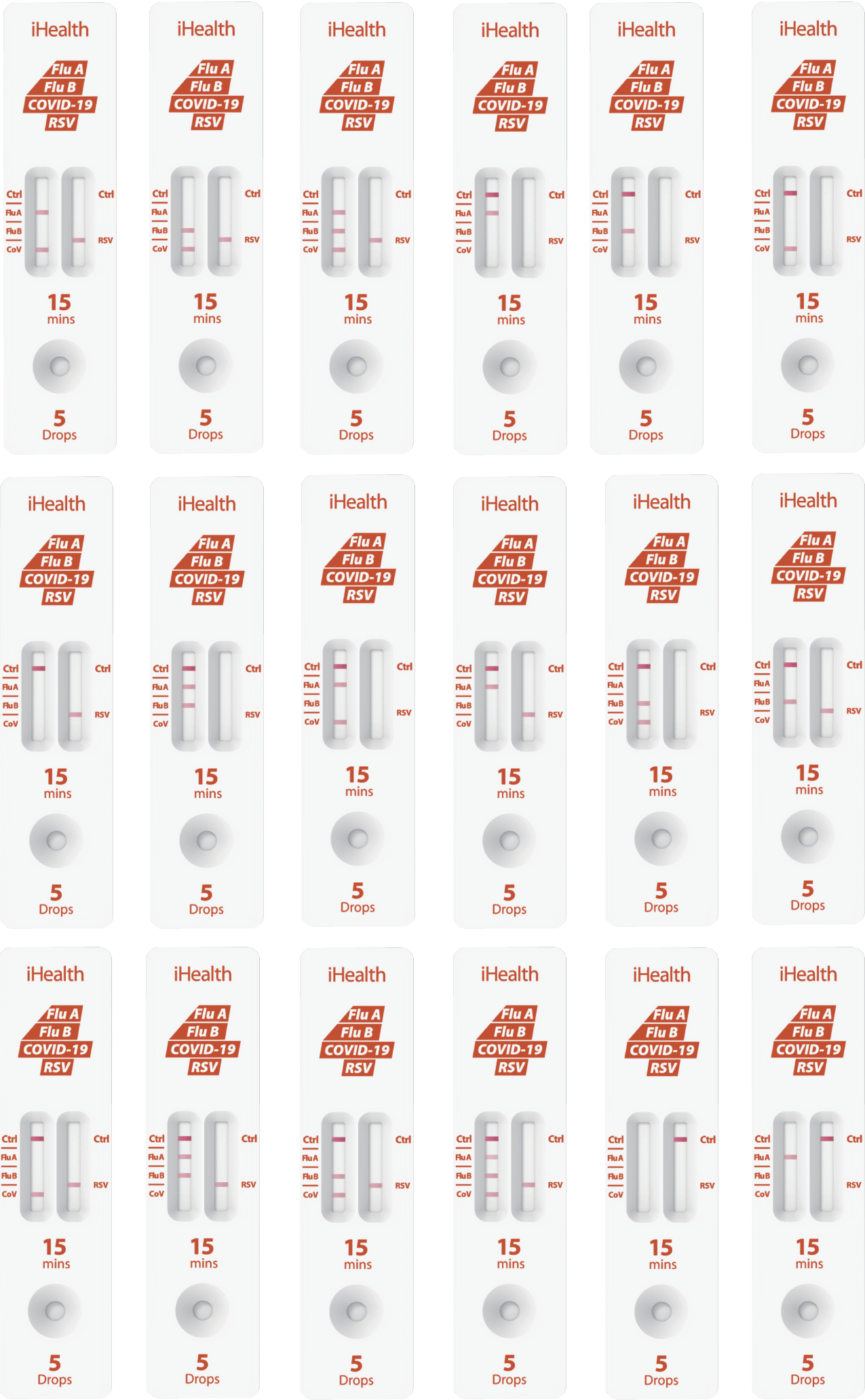


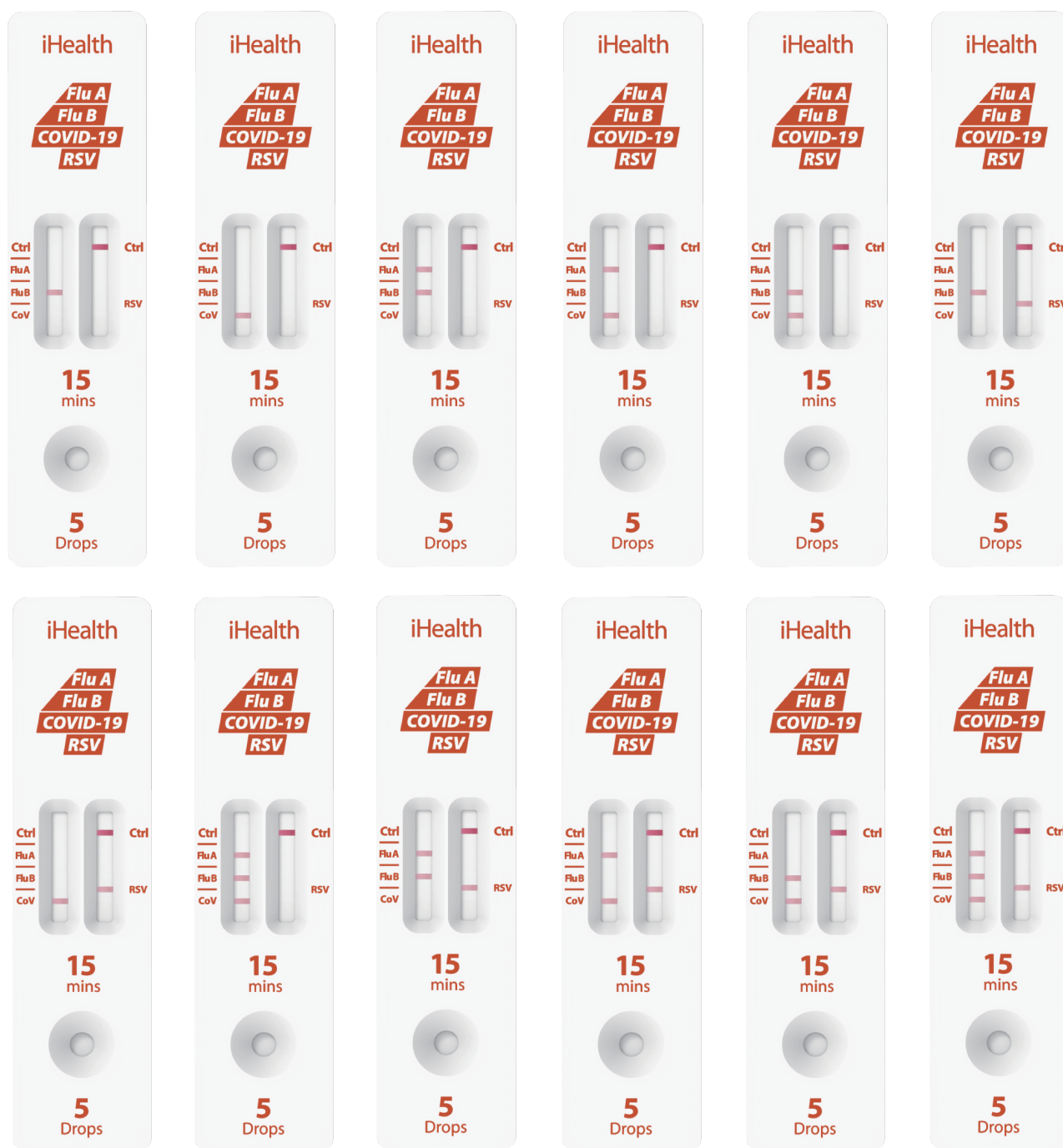


Missing 'Ctrl' line on ONE or BOTH strips


Invalid Result: Check to see if a line is visible at the control line 'Ctrl' on both strips. If one or both 'Ctrl' lines are missing, the result is **INVALID**, do not continue reading the results. Repeat the test with a new sample and a new test kit.

Some causes of invalid results are because of not following the directions correctly or the test may have deteriorated beyond the expiration date. It is recommended a specimen be re-tested using a new test.

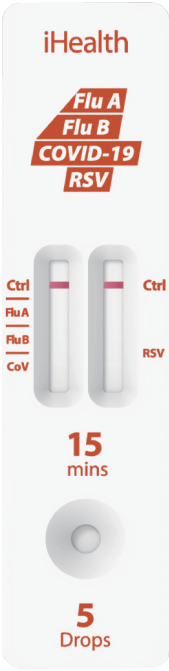


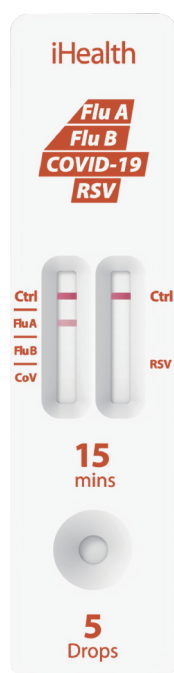


Flu A

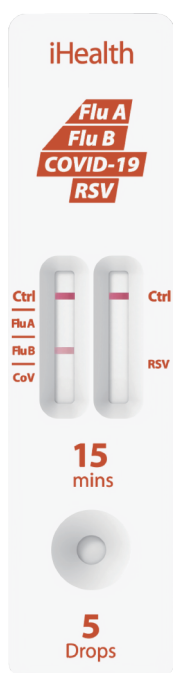

Flu B

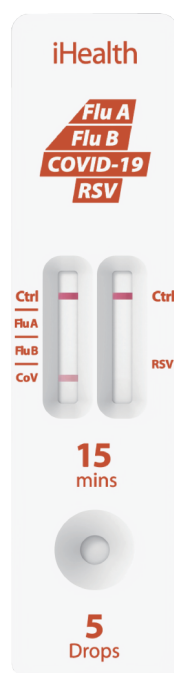


COVID-19

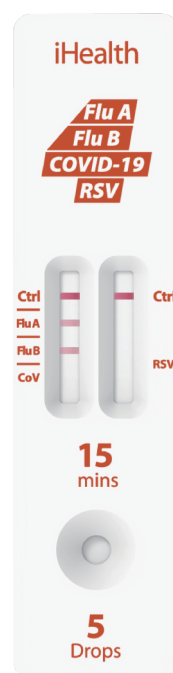


Flu A & Flu B




Flu A & COVID-19

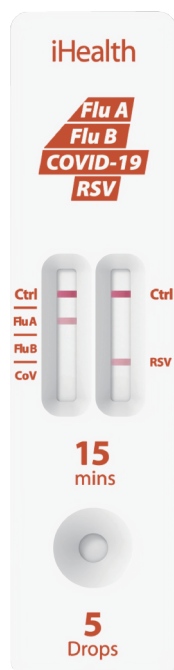


Flu A & RSV

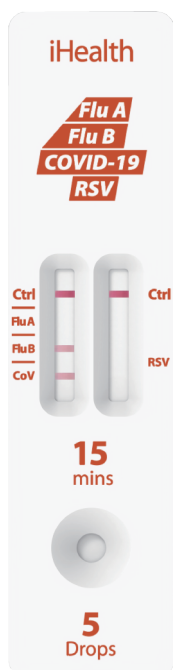


Flu B & COVID-19

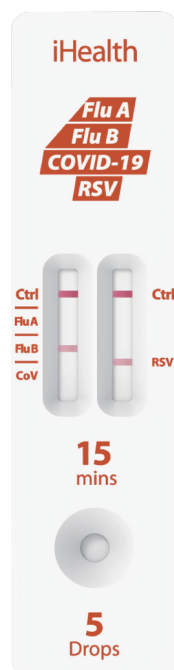


Flu B & RSV

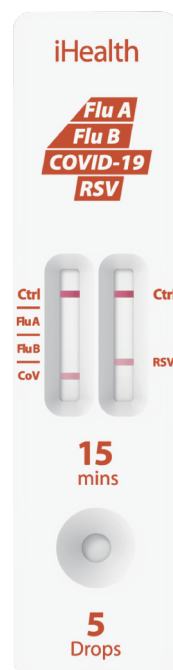


COVID-19 & RSV

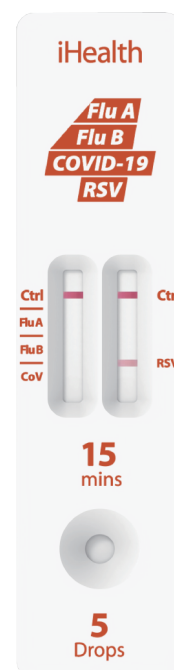


RSV

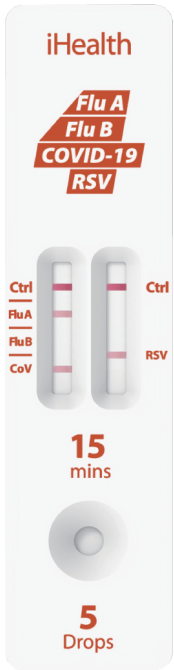


Flu A & COVID-19 & RSV

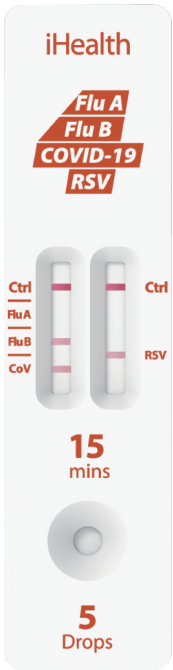


Flu B & COVID-19 & RSV

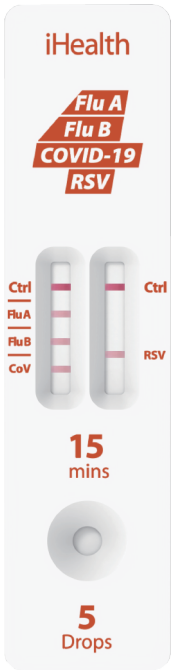

Flu A & Flu B & COVID-19 & RSV

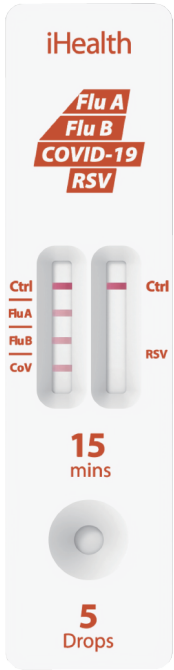

Flu A & Flu B & COVID-19

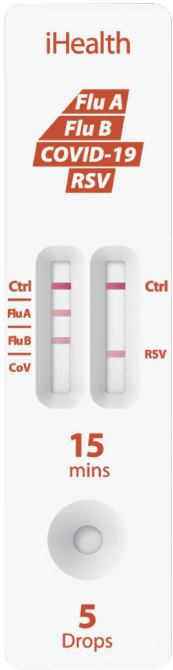

Flu A & Flu B & RSV